

On·Call® MultiPro System

Hemoglobin A1c Test



Operation Steps

Preparation

Switch on the Analyzer.

Refer to the following user instructions for detailed information:

- On Call® MultiPro System Instruction Manual
- On Call® MultiPro HbA1c Test Kit Package Insert
- On Call® MultiPro HbA1c Controls Package Insert

For a Patient Test

Sample Collection




Or



Fill the Capillary completely with blood from finger.

Fill the Capillary completely with blood in blood collection tube.


Touch only the tip of the Capillary to the blood sample on the stopper.

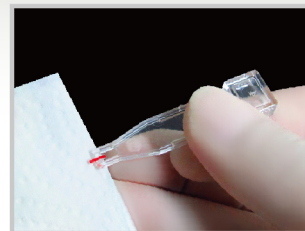


Use a lint-free tissue to carefully wipe the outside of the Capillary.

Note: To avoid excess blood on the Sampling Device, DO NOT submerge the Sampling Device in the blood collection tube. Inspect the Capillary for the presence of bubbles. If bubbles are present, discard the Sampling Device and repeat the procedure using a new Sampling Device.

Performing a Test

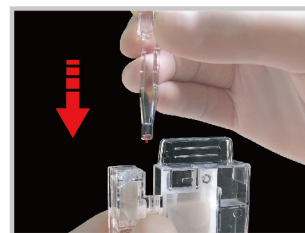


Enter Patient ID/Patient Name/Sample ID before the test. Snap the filled Sampling Device into Test Cartridge.

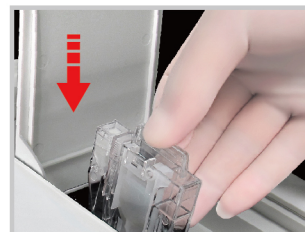

Open the cartridge compartment cover. Insert the Test Cartridge into the cartridge compartment with the label facing left.



Close the cover to start the test. Results will be displayed on the screen after 5 minutes.

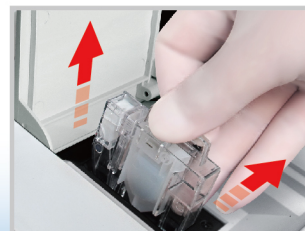

Open the cartridge compartment cover. Remove the used Test Cartridge. Close the cartridge compartment cover.

Note: Gently push the tab on the cartridge to the right when remove the cartridge.



Specification

Feature	Specification
Feature	Specification
Methodology	Immunoturbidimetry
Specimen Volume	1µL
Operation Method	One step operation
Test Time	About 6 minutes/test
Application	Only use <i>On Call® Multipro</i> reagents
Measurement Range	3.8% to 15.0% (17.8mmol/mol-141mmol/mol)
Specimen	Capillary (fingertip) and venous whole blood
Units of Measure	NGSP (%) & IFCC (mmol/mol)
Memory	5000 test results
Automatic Sleep Mode	30 minutes
Analyzer Size	157 x 313 x 158 mm
Display Size	TFT 1024 x 600 touch screen
Weight	2.3kg without power supply and printing paper
Analyzer Storage Conditions	-20 - 50°C; 10%-90%RH
Test Cartridges Storage Temperature & Shelf Life	2 - 8°C (35.6 – 46.4°F) Shelf Life: 24 months
	9 - 30°C (48.2 - 86°F) Shelf Life: 3 months
Operating Conditions	15 - 30°C (59 - 86°F) ; 10-90%RH
Printer	Inner thermal printer; External printer interface
Data Transport	RS232, USB, Internet
Power Source	100 – 240 VAC

Catalog

Catalog No.	Product Name	Contents
G116-112	<i>On Call® MultiPro</i> Analyzer	1 Analyzer Printer Paper Rolls Power Cord Power Supply Adapter Quick Reference Guide Instruction Manual
G116-113	<i>On Call® MultiPro</i> Analyzer (with Barcode Reader)	1 Analyzer Printer Paper Rolls Power Cord Power Supply Adapter Quick Reference Guide Instruction Manual Barcode Reader
G136-112	<i>On Call® MultiPro</i> HbA1c Test Kit	10 Test Cartridge 1 Sampling Device Package Insert
G126-122	<i>On Call® MultiPro</i> HbA1c Controls	5 x 1µg HbA1c Control Level 1 5 x 1µg HbA1c Control Level 2 1 HbA1c Control ID Cartridge Package Insert
G126-123	<i>On Call® MultiPro</i> Printer Paper Rolls	
G126-124	<i>On Call® MultiPro</i> Barcode Reader	
G126-125	<i>On Call® MultiPro</i> Data Transfer Cable	